

Oxygen (%sat)

Fig 4.3. Water quality index (Note- If dissolved oxygen is greater than 140%, the quality index equals 50)

### BIOCHEMICAL OXYGEN DEMAND (BOD)

(The Measure of Water Pollution)

To measure the purity of water BOD (Biochemical oxygen demand) unit is used. **BOD is the quantity of oxygen which is required by micro-organisms to decompose harmful substances present in one liter water.** BOD value of pure drinking water must be 1 ppm. When oxygen is not available according to BOD micro-organisms fail to decompose harmful substances, then these substances are decomposed by anaerobic method and again some other harmful substances are formed, as a result of this water gives bad smell..

The BOD test is the most important measure of the polluting capacity of organic effluents. BOD is measured by keeping a sample of water containing known amount of oxygen for five days at 20°C in the dark. At the end of this period the oxygen content is again measured. A high BOD indicates intense level of microbial pollution.

BOD measures the potential of the organic matter in a sample to deplete the avail-



able dissolved oxygen. The degradation of organic matter by the microorganisms present in a water body requires oxygen. This can quickly deplete the available (Dissolve) oxygen in the water. When the dissolved oxygen levels drop too low, many aquatic species perish. In fact, if the oxygen level drops to zero, the water will become septic. When organic compounds decompose without oxygen, it gives rise to the undesirable odours usually associated with septic or putrid conditions.

Biological Oxygen Demand (BOD) is one of the most common measures of pollutant organic material in water. BOD indicates the amount of putrescible organic matter present in water. Therefore, a **low BOD is an indicator of good quality water**, while a **high BOD indicates polluted water**.

DO is the actual amount of oxygen available in dissolved form in the water. When the DO drops below a certain level, the life forms in that water are unable to continue at a normal rate. The decrease in the oxygen supply in the water has a negative effect on the fish and other aquatic life. Fish kills and an invasion and growth of certain types of weeds can cause dramatic changes in a stream or other body of water. Energy is derived from the oxidation process.

Table 4.4: BOD level of different types of water.

BOD Level (in ppm)	Water Quality
1 - 2	<b>Very Good</b> There will not be much organic waste present in the water supply.
3 - 5	<b>Fair - Moderately Clean</b>
6 - 9	<b>Poor - Somewhat Polluted</b> Usually indicates organic matter is present and bacteria are decomposing this waste.
100 or greater	<b>Very Poor - Very Polluted</b> Contains organic waste.

{Note: Generally, when BOD levels are high, there is a decline in DO levels. This is because the demand for oxygen by the bacteria is high and they are taking that oxygen from the oxygen dissolved in the water. If there is no organic waste present in the water, there won't be as many bacteria present to decompose it and thus the BOD will tend to be lower and the DO level will tend to be higher. At high BOD levels, organisms such as macroinvertebrates that are more tolerant of lower dissolved oxygen (i.e. leeches and sludge worms) may appear and become numerous. Organisms that need higher oxygen levels (i.e. caddisfly larvae and mayfly nymphs) will NOT survive}.

BOD specifies the strength of sewage. In sewage treatment, to say that the BOD procedure that measures the amount of oxygen consumed by living organisms while they are utilizing the organic matter present in waste, under conditions similar in nature.

#### Measurement of B.O.D.

The BOD test takes **5 days to complete** and is performed using a dissolved oxygen test kit. The BOD level is determined by comparing the DO level of a water sample taken immediately with the DO level of a water sample that has been incubated in a dark oxygen required for the decomposition of any organic material in the sample and is a good approximation of the BOD level.

1. Take 2 samples of water
2. Record the DO level (ppm) of one immediately using the method described in the dis-



- solved oxygen test.
- Place the second water sample in an incubator in complete darkness at 20 °C for 5 days. If you don't have an incubator, wrap the water sample bottle in aluminum foil or black electrical tape and store in a dark place at room temperature (20 °C or 68 °F).
  - After 5 days, take another dissolved oxygen reading (ppm) using the dissolved oxygen test kit.
  - Subtract the Day 5 reading from the Day 1 reading to determine the BOD level. Record your final BOD result in ppm.


### CHEMICAL OXYGEN DEMAND (COD)

In environmental chemistry, the **chemical oxygen demand (COD)** test is commonly used to indirectly measure the amount of organic compounds in water. Most applications of COD determine the amount of organic pollutants found in surface water (e.g. lakes and rivers), making COD a useful measure of water quality. It is expressed in milligrams per liter (mg/L), which indicates the mass of oxygen consumed per liter of solution.

The basis for the COD test is that nearly all organic compounds can be fully oxidized to carbon dioxide with a strong oxidizing agent under acidic conditions. The amount of oxygen required to oxidize an organic compound to carbon dioxide, ammonia, and water is given by-




This expression does not include the oxygen demand caused by the oxidation of ammonia into nitrate. The process of ammonia being converted into nitrate is referred to as *nitrification*. The following is the correct equation for the oxidation of ammonia into nitrate.

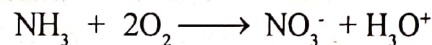

The second equation should be applied after the first one to include oxidation due to nitrification if the oxygen demand from nitrification must be known. Dichromate does not oxidize ammonia into nitrate, so this nitrification can be safely ignored in the standard chemical oxygen demand test.

### Measurement of COD

The test of COD is carried out to measure the contents of organic matter of sewage and natural water. Most types of organic matter are oxidized by a boiling mixture of chromic and sulfuric acids. A sample is refluxed (Boiled) for two hours in strongly acid solution with a known excess of potassium dichromate ( $K_2Cr_2O_7$ ) in presence of a catalyst  $AgSO_4$  and  $HgSO_4$ . It is then cooled and the remaining unreacted  $K_2Cr_2O_7$  is titrated with ferrous ammonium sulphate to determine the amount of  $K_2Cr_2O_7$  consumed and the oxidizable matter is calculated in terms of oxygen equivalent.

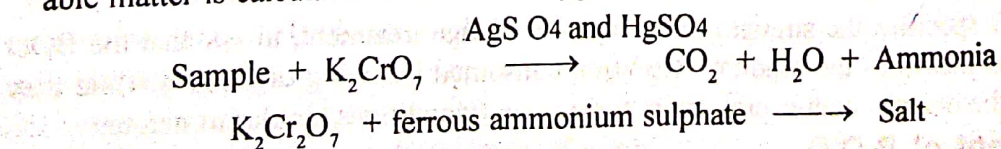

### Calculations

The COD is calculated with the help of following formula

$$COD \text{ mg/l} = \frac{(V_1 - V_2) \times 8 \times 1000}{x}$$

Where :  $V_1$  = Initial volume of ferrous ammonium sulphate  
 $V_2$  = Final volume of ferrous ammonium sulphate



$X$  = Volume of sample taken.

### **Advantage of COD over BOD**

COD is usually carried out to determine the pollutional strength of sewage because of its several advantages over BOD test. These advantages are-

1. It takes relatively less time (Only 2 hours).
2. Industrial waste do not responds to BOD test, so COD test is a must.
3. If toxic materials are present in sewage, they are likely to interfere with BOD, in such cases COD is found to be very useful.